

### **POSTER PRESENTATION**

Open Access

# Preliminary studies of cell culture strategies for bioprocess development based on HEK293 cells

Leticia Liste-Calleja\*, Jonatan López-Repullo, Martí Lecina, Jordi Joan Cairó

From 23rd European Society for Animal Cell Technology (ESACT) Meeting: Better Cells for Better Health Lille, France. 23-26 June 2013

#### **Background**

The use of human embryonic kidney cells (HEK293) for recombinant protein or virus production has gained relevance along the last years. They are specially recommended for transient gene expression and adenovirus or adeno-associated virus generation [1,2]. To achieve high volumetric productivities towards bioprocess optimization, the concentration of biocatalizer (i.e. animal cells) must be enhanced. The limits for cell growth are mainly related to the accumulation of metabolic by-products, or the depletion of nutrients [3]; therefore, cell cultures strategies must be developed. In this work, we have explored Punctual Feeding and Media Replacement cell culture strategies to over perform the limit on Xvmax encountered on batch culture mode. Finally, we scaled up cell culture in order to control other parameters (i.e. pO<sub>2</sub>) that could be limiting cell growth.

#### Materials and methods

The cell line used in this study was HEK293SF-3F6 (kindly provided by Dr. A.Kamen, NRC-BRI). The basal medium for all cell cultures was SFMTransFx-293 (Hyclone, Thermo Scientific) supplemented with 5% (v/v) of FBS and 4 mM GlutaMAX (Gibco, Invitrogen). For Punctual Feeding and FedBatch Fementation Cell Boost 5 (Hyclone, Thermo Scientific) was used. Batch, media replacement and punctual feeding experiments were performed in 125-ml plastic shake flasks (Corning Inc.) shaken at 110 rpm in an orbital shaker at 37°C, 95% humidity, 5%  $\rm CO_2$  incubator. FedBatch Fermentation was carried out in Bioreactor Braun-MCD (2 L) with mechanical agitation at 80 rpm, pH set point 7.1 and pO<sub>2</sub> set point 50%. Viable cell density and viability were determined by trypan blue exclusion method and


manual counting using a haemocytometer. Glucose and lactate were analysed in an automatic analyser, YSI (Yellow Springs Instrument, 2700 Select).

#### **Results**

Characterization of HEK293 cell culture in batch operation was initially performed. It was encountered that cell growth was extended for 168 h reaching approximately  $7.10^6$  cell/mL of cell density (Figure 1.1). Nevertheless, maximal cell growth rate  $(\mu_{max})$  was only maintained for 96 h. As glucose and lactate were not at limiting concentrations [4], nutrient limitation different from glucose arose as the first hypothesis for this decrease on cell growth rate. Therefore, punctual additions of nutritional supplement for HEK293 were carried out. Xv<sub>max</sub> was significantly increased in comparison to basal media, reaching cell densities as high as  $17 \cdot 10^6$  cell/mL (Figure 1.2). Nevertheless, we could not overcome this limit on  $Xv_{max}$ regardless the number of punctual feedings performed. Moreover, nutrient addition did not elongate  $\mu_{max}$  period (t<sub>11</sub>). These results suggested that by-product accumulation different from lactate could be limiting cell growth. In order to validate the hypothesis, complete media replacement (up to three replacements) was studied (Figure 1.3). Although this strategy enabled to extend t<sub>u</sub> up to 168 h of cell culture, the maximal cell density reached was similar to nutrient addition strategy (1MR: 12·10<sup>6</sup> cell/mL; 2MR: 16·10<sup>6</sup> cell/mL; 3MR: 18·10<sup>6</sup> cell/ mL). This limit on Xv<sub>max</sub> encountered on shake flask might be related to a limitation on pO2. Thus, the cell culture system was changed towards a bioreactor with controlled pO2 (maintained between 20-60% of air saturation). In addition, a continuous feeding using a pre-fixed profile addition was implemented. As it can be noticed in Figure 1.4, FedBatch operation in bioreactor enabled to beat the limit encountered in shake flask system, reaching cell densities of 27·10<sup>6</sup> cell/mL.

<sup>\*</sup> Correspondence: Leticia.Liste@uab.cat Chemical Engineering Department, Universitat Autònoma de Barcelona, Cerdanyola del Vallès, 08193, Spain



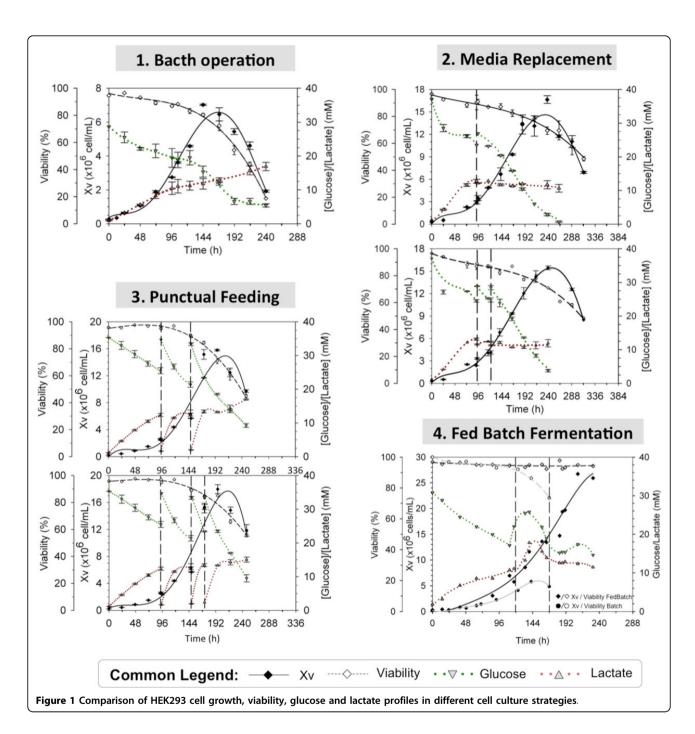

#### **Conclusions**

Punctual feeding and media replacement overcame the limit of  $7{\cdot}10^6$  cell/mL encountered in batch mode operation indicating that nutrient depletion was one of the causes of that limit. Nevertheless, the elongation of  $t_{\mu}$  found out performing MR suggests that the accumulation of by-products might not be ruled out.

The new limit on  $Xv_{max}$  ( $\approx 17-18\cdot 10^6$  cell/mL) encountered regardless the cell culture strategy, was outperformed by transferring  $O_2$  more efficiently in

bioreactor system, reaching cell densities as high as  $Xv_{max} = 27 \cdot 10^6$  cell/mL. The monitoring and control of cell culture parameters (i.e.  $pO_2$ , pH) will enable to develop more accurate feeding strategies in order to achieve higher cell densities than those presented here (on going work).

#### Acknowledgements

We would like to thank Dr. Amine Kamen (BRI-NRC, Canada) for kindly providing the HEK 293 cell line.

#### Published: 4 December 2013

#### References

- Nadeau I, Kamen A: Production of adenovirus vector for gene therapy. Biotechnology advances 2003, 20:475-489.
- Geisse S, Fux C: Recombinant protein production by transient gene transfer into Mammalian cells. Methods in Enzymology 2009, 463:223-238.
- Butler M: Animal cell cultures:recent achievements and perspectives in the production of biopharmaceuticals. Appl Microbiol Biotechnol 2005, 68:283-291.
- Petiot E, Jacob D, Lanthier S, Lohr V, Ansorge S, Kamen A: Metabolic and Kinetic analyses of influenza production in perfusion HEK293 cell culture. BMC Biotechnol 2011, 11:84-96.

#### doi:10.1186/1753-6561-7-S6-P52

Cite this article as: Liste-Calleja *et al.*: Preliminary studies of cell culture strategies for bioprocess development based on HEK293 cells. *BMC Proceedings* 2013 **7**(Suppl 6):P52.

## Submit your next manuscript to BioMed Central and take full advantage of:

- Convenient online submission
- Thorough peer review
- No space constraints or color figure charges
- Immediate publication on acceptance
- Inclusion in PubMed, CAS, Scopus and Google Scholar
- Research which is freely available for redistribution

Submit your manuscript at www.biomedcentral.com/submit

